PLUS⁺

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declares that the object of the declaration is the product: Massage device

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conforms to essential requirement according to ECC Directive EMC 2014/30/UE, LVD 2014/35/UE, ROHS 2011/65/UE, ERP Regulation 1275/2008 and subsequent modifications and DECRETO 15 October 2015, n. 206 applicable to this product. In addition, the object of the declaration described above, is in conformity with the relevant Union harmonization legislation:

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This declaration of conformity is issued under the sole responsibility of the manufacturer.

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Additional information:

Torino, 23/09/2019

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Signed for and on behalf of:

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ENGLISH

Dear customer,

thank you for choosing the device Mesis® PressoMassaggio® Plus+. You'll appreciate the high quality of the materials used as well as the features that set this model apart.

Along with wishing you pleasant use of the product, we kindly remind you to read the manual, the contraindication and the warnings before use.

The product must be used according to what is written in the present manual.

The producer and the distributor shall not assume any liability for any different use other than the one indicated in the present manual.

Laws and regulations, to which reference is made in this manual, are European and Italian.

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Via Pigafetta 65, 10129 Torino

fax (+39) 011.0432817

email: info@mesis.eu

Reproduction, even in part, of any photo or text contained within this manual is prohibited.

Index

Warnings.....	30
Features of Mesis® PressoMassaggio® Plus+	38
Mesis® PressoMassaggio® Plus+ Mod. PSM1000P	40
Mesis® PressoMassaggio® Plus+ operating instructions.....	42
Applicators and accessories.....	44
Instructions for use applicators with 4 chambers.....	46
Maintenance and troubleshooting.....	49
Warranty	50
EU declaration of conformity	51

Warnings

Before using the device, in case of doubts, consult your doctor and read the instructions attentively for use and, in particular, safety warnings, thereby adhering to these.

Conserve the present manual for the entire service life for reference.

In the case of transferring the product to third parties, please also provide the documentation in its entirety.

After having removed the device from the packaging, check to ensure that the scope of supply is complete and check for any possible transport damage. If in doubt, do not use the device and contact the authorised service centre.

The packaging material is not a toy for children! Keep the plastic bags away from children; danger of suffocation!

Before connecting the device, check that the voltage data shown on the identification plate correspond to the specifications of electric mains available.

The identification label is located on the device.

The present device must only be used for the intended purpose it was designed for, or as a massage unit. Any other use is considered improper and therefore dangerous.

In order to ensure the best level of protection, it is advised to install a differential current device the rated operating differential current of which does not exceed 30mA (ideally 10mA). Ask your installer for advice concerning this.

Caution:

Do not use the present device near bathtubs, sinks or other receptacles containing water!
Do not submerge the device into water!

The appliance can't be used by children aged 8 or less and by persons with reduced physical, sensory or mental abilities, or lacking experience or necessary knowledge, provided they are under surveillance or after they have received instructions regarding the safe use of the appliance and the understanding of the dangers inherent to it. Children must not play with the appliance. Cleaning and maintenance intended to be carried out by the user must not be carried out by children without supervision. Keep the appliance and its cord out of reach of children under 8 years.

DO NOT use the device with wet hands, damp or bare feet.

DO NOT pull the power cable or the device itself when disconnecting it from the power outlet.

DO NOT expose the device to moisture or the influences of weather (rain, sun).

Do not block the air openings of the appliance and do not place the appliance on a soft surface such as a bed or couch, where the openings may become blocked.

Keep the air openings free of any obstructions.

In the case of breakdown or malfunction of the device, turn it off and do not tamper with it.

Only contact an authorised technical service centre for any possible repairs.

If the power cord is damaged, it must be replaced by an authorised technical service centre in order to prevent any risk.

Use a moist cloth to clean the unit and check the contents of the present manual.

Please note that the apparatus must only be used indoors and in accordance with the operating instructions. If the power cord is faulty, it must be replaced by a VIMA S.r.l. Società a Socio Unico customer service centre.

The electrical unit to which the device is connected has to be in conformity with the Ministerial Decree No. 37/08 (Italian Rules, then check the applicable regulations in your country) and, in particular, all the safety measures have to be in place for the unit - thermo-magnetic switches on the main power lines.

Please store and handle the product with care. It is prohibited to use the product if damage is present.

Do not position the device or its accessories near sources of heat.

Only use in domestic environments.

Photos and features the product refer to what is currently printed in the manual, differences of minor importance may be possible during the course of production. Possible differences shall have no influence on the safety or on the performance of the product.

In order to avoid flame or shocks of various nature, do not expose the device to rain or liquid of any type.

The manufacturer shall provide information on the correct use of the items. Therefore, always adhere to these instructions that assure the correct and appropriate installation as well as functionality. For this reason, always keep these instructions for future reference.

Position the switch to 0 (zero) or OFF after discontinuing use.

All of the device's technical specifications must be borne in mind. Check the symbols that are on the ID tag of each item and the icons visible in the section (page 33-34-35-36-37) of the safety instructions.

This appliance may only be used indoors.

It is not suitable for installation in bathrooms.

The item is only suited to be directly connected to the mains using a proper power supply.

Electrical socket: in the case of damage, the electric cable must only be replaced by the manufacturer, the distributor or by an expert technician in order to prevent risks.

During installation ensure you do not make holes near the existing electrical wiring in the wall or ceiling.

Connecting cables must never be twisted and squashed into the space between the device and the supporting surface.

This product isn't suitable for children.

Do not use in the following cases:

- subjects with capillary fragility, arterial hypertension or heart disease;
- subjects with a pacemaker or an internal defibrillator.
- carriers of active implants such as pacemakers, insulin pumps or other devices;
- carriers of metallic prostheses;
- subjects with inflammatory processes in place, with skin lesions, with neoplasms;
- pregnant women;
- subjects with sensitive skin, who may be restricted to powers of very low disbursement;
- apply only on undamaged skin;
- never use the product directly on swollen or inflamed areas or skin rashes.

In case of doubts, consult with your medical specialist before use.

How to operate and further warnings:

The application time varies depending on how the massage is carried out and, as a rule, this can be vary between 15 and 30 minutes.

Other than those listed above, attentively follow the indications, warnings and cautions for use specified within the present manual.

Mesis® PressoMassaggio® Plus+ machine is designed, manufactured and tested in compliance with the most stringent European safety standard regulations.

Any damage caused by using the device in exceptional circumstances (use outdoors, in wet areas, wet zones and exposure to direct sunlight ...) is not covered by the warranty terms.

The period of warranty shall begin from the date of purchase and shall only be deemed valid if a document of purchase, invoice or receipt is presented. The warranty shall lapse if the item has not been installed according to the specified instructions or if it is repaired or tampered with.

Furthermore, the manufacturer shall not assume any liability for damages caused due to misuse or incorrect application of the instructions.

We shall reserve the right to perform technical modifications or modifications to the design.

Environmental warnings

In compliance to Legislative Decree No. 49 of March 14, 2014 (Italian regulatory reference).

„Implementation of Directive 2012/19/EU on waste electrical and electronic equipment (WEEE)“

The crossed bin symbol on the appliance indicates that the product at the end of its useful life must be collected separately from other waste. The user must therefore confer the equipment complete with essential components at the end of their life to the appropriate collection centers for electronic and electrotechnical waste, or return it to the seller when new equipment of an equivalent type is purchased, in one by one, or one zero for equipment with major side lower than 25cm. The adequate separate collection for the subsequent start of the equipment discharged to recycling, massage and compatible environmental disposal helps to avoid possible negative effects on the environment and health and favors the recycling of the materials that make up the equipment. The illegal disposal of the product by the user involves the application of administrative sanctions in compliance to Legislative Decree No. 49 of March 14, 2014 (Italian regulatory reference). It is advisable to check the regulations in force in your country regarding the disposal of waste and Environmental warnings.

Cleaning

Always disconnect the electric power plug before cleaning or maintenance and in the case of non-use of the device.

Do not use stringent or abrasive detergents for cleaning the surface: use a damp cloth for cleaning the device when necessary.

If dirt is difficult to remove, use a light detergent and delicately rub the area in question.

The parts that come into contact with the skin must be clean between massages etc. using normal cleaning methods: gentle disinfectant on a soft cloth.

Disconnect from the power supply before cleaning.

It is recommended to properly sterilise and/or disinfect all the parts that come into contact with the individual undergoing.

Clean items with a damp cloths and do not use solvents or abrasive products.

To clean the appliance, use a damp cloth and check the contents of this manual.

Avoid spreading liquids on all electrical parts.

Warnings

- Prior to use, read each section of the manual attentively;
- Store this manual close to the device so that you can refer to it at any time;
- This manual explains the correct use of the devices in detail in order to prevent dangers, personal harm and material damage.

Meaning of the symbols specified in the present manual

 Attention	 Precautions
<p>If these instructions are not adhered to, this could cause severe and fatal injuries.</p>	<p>If these instructions are not adhered to, the user could cause injuries or damage to the device.</p>

Indications

 It indicates a restriction.	 It is an indication to be correctly followed.
 Disassembly is prohibited.	 Indicates to disconnect the power cable of the wall outlet.

Attention

Contraindications

- Individuals that suffer from inflammation or injuries that are at risk of bleeding, both internally or externally may not use the product including those suffering from the discharge of pus.
- In the case of general discomfort, flu or illness, do not use the product.
- Pulmonary oedema
- Thrombophlebitis
- Heart failure
- Infections
- Deep vein thrombosis
- Pulmonary embolism episodes
- Injuries, lesions or tumours at the area of application or nearby where venous and lymphatic return could result in injury
- carriers of active implants such as: pacemakers, insulin pumps or other devices;
- Lymphangitis
- Plexopathy and neuropathy
- Carriers of metallic prostheses
- Subjects with inflammatory processes in place, with skin lesions, with neoplasms
- Pregnant women

- Subjects with sensitive skin, who may be restricted to powers of very low disbursement
- apply only on undamaged skin.

Consult your doctor to check whether you have any disease which is incompatible with the Mesis® PressoMassaggio® Plus+ massage. A non-exhaustive list of incompatible diseases are listed in the „Contraindications“ section.

Warnings

Attention

! **If you feel bad during use or you have any pain, do not use the device.** Immediately suspend its use.

! **For the first time, use the device for less than 20 minutes.** Possibly increase the time of the following sessions based upon your body's own tolerability. Start the massage at a minimum level of pressure.

⊘ **It is prohibited for children to use the device.** Keep out of the reach of children. Danger of electric shock or other damage.

⊘ **Do not use water or other substances that could seep into the inside of the device.**

⊘ **In the event of breakdown or malfunction, do not open or try to repair the product.** Possible electric shocks, personal injury and material damage.

⊘ **Your hands must be dry when connecting or disconnecting the power cable from the electric mains (outlet).** Possible electric shocks, personal injury and material damage. Fire danger.

! **Use only 220V single-phase current. After using the device, properly disconnect the cable from the electric mains.** Possible electric shocks, personal injury or material damage. Fire danger.

⊘ **Keep the power cable free of any obstructions; don't place anything on top of it, crush it or pinch it with any object.** Possible shock or electrical leakage. Fire danger.

Warnings

Attention


Use an electrical outlet exclusively for the device. Too many plugs can cause a fire, electrical leakage or electric shocks.


Do not hit the device; do not let it fall. This could cause serious damage to the device.

Do not insert any object or product into the device's interior. Risks of damage to the device and personal injury.

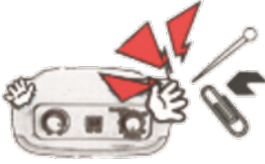Do not place the device or its accessories near sources of heat. These could cause serious damage to the device and/or its accessories.

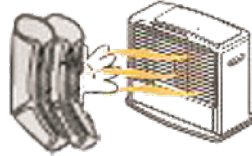When wearing the legs cuff, the abdominal belt or arm cuff, do not move or walk.

Do not subject the device and its accessories to direct sunlight or gas heaters. They might cause decay, deformation and discolouration of the products.

Do not use the device for any other use than massage. Risks of personal injury or damage to the device.

Do not bend, knot, cut or break the electric cable.

Warnings


Attention

- Do not place any object to the top of the device.** Do not put weight on the unit or on its connections, nor on electric cables and connecting air pipes. This could cause variations in pressure.

- If the device has not been used for a longer period of time or is been exposed to long-term cold,** it may not function regularly. Therefore, once it is turned on, allow it to operate at a low level of pressure for at least three inflation cycles. Do not leave the unit in a humid place or in the rain.

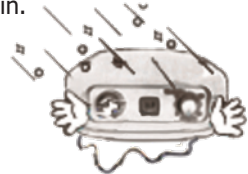- To clean the accessories (legs cuff, arm cuff etc.) do not use water, oil, gasoline or chemical substances. Use a damp cloth.** This could ruin it, thereby shortening its service life.

- Do not use the accessories directly on the skin.** Under the applicators, wear appropriate underwear.

- To clean the device, do not use water, oil, gasoline or chemical substances. Use a damp cloth.**

- Do not bend, knot or cut the air tubes.**

Features of Mesis® PressoMassaggio® Plus⁺

Main features

- Mesis® PressoMassaggio® Plus⁺ is sequential and intermittent;
- Adjustable continuous pressure
- Optimal massage by means of the 4 cuffs;
- Integrated timer to adjust the duration of massage;
- Control valve can supply a constant pressure;
- Fast connection of accessories and easy use of the unit.


Terminology

- **Program:** method of inflating the cuffs;
- **Chambers or segments:** elements that the applicator is made up of and that fill with air. In the Mesis® PressoMassaggio® Plus⁺ the cuffs of the applicators are 4.
- **Cycle:** one cycle equals the inflation of all cuffs, from the first to the last one independent of the program used. After inflation of the last cuff, there is a break and then a new cycle starts;
- **Applicator:** pneumo-soothing elements composed of several cuffs that are sequentially inflated and regulated by the selected program, (legs cuff, arm cuff and abdominal/buttocks belt);
- **Time:** duration of session (adjustable from 1 to 30 minutes);
- **Pressure:** quantity of air introduced into the chambers, the pressure is adjustable and continuous;
- **Connector:** it is an accessory made up of joined rubber tubes and allows the applicators to be connected to the machine or to the splitter;
- **Splitter:** it is an accessory made of rubber tubes joined together. The use of 1 splitter, allows the use of the two leg cuffs while the use of 2 splitters, allows the use of three accessories at the same time, for example the two leg cuffs and the abdominal/buttocks belt;
- **Nozzle:** end part (air outlet) of each connector's hose that must be inserted into the seats located on the applicators;
- **Seat:** located on the applicator. This is where the connector's nozzles must be inserted.

Program (sequence)

Program A

The pressure is exerted by one cuff at a time.

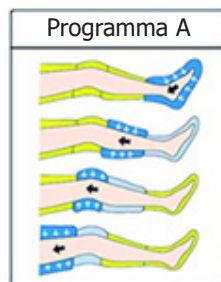Inflation direction

Arm cuff:
inflation direction

Abdominal/buttocks belt:
inflation direction

Legs cuff:
inflation direction

In the above images, the arrows indicate the correct inflation of the applicators:

- if the arm cuff is used, the inflation starts with the hand to end the shoulder;
- If the abdominal / buttocks waist is used, the inflation begins at the top of the abdomen and proceeds downwards;
- if you use your legs cuff, the inflation begins with the foot to end the groin.

Mesis® PressoMassaggio® Plus⁺ Mod. PSM1000P

Main technical features

- 1 program
- Applicators with 4 chambers
- Pressure regulation 30~240 mmHg \pm 10%
- Adjustable massage time 0~30 min.
- Power supply: 220 – 240 Vac, 50/60Hz
- Electrical consumption: 25W
- Dimensions: 220X200X115mm
- Device weight: 2 kg
- Warranty: 24 months
- Declaration of conformity CE

Configurations

Cod. PSM1000P-1B	Cod. PSM1000P-F	Cod. PSM1000P-1G
		
<p><i>Applicators:</i> 1 arm cuff <i>Accessories:</i> 1 single connector</p>	<p><i>Applicators:</i> 1 belt <i>Accessories:</i> 1 single connector</p>	<p><i>Applicators:</i> 1 leg cuff <i>Accessories:</i> 1 single connector</p>
Cod. PSM1000P-2G	Cod. PSM1000P-2GK	Cod. PSM1000P-2GKB
		
<p><i>Applicators:</i> 2 leg cuffs <i>Accessories:</i> 2 single connectors and 1 splitter</p>	<p><i>Applicators:</i> 2 leg cuffs and 1 belt <i>Accessories:</i> 3 single connectors and 2 splitters</p>	<p><i>Applicators:</i> 2 leg cuffs, 1 belt and 1 arm cuff <i>Accessories:</i> 4 single connectors and 2 splitters</p>

Mesis® PressoMassaggio® Plus⁺ operating instructions

IMPORTANT: to prevent the compressor (heart of the equipment) from being subjected to sudden movements and blows during transport, it has been blocked. This solution was adopted to guarantee the integrity, functionality and quality of the equipment. The compressor must then be unlocked in the following mode.


BEFORE CONNECTING THE POWER CORD AND PERFORMING ANY OPERATION ON THE MACHINE, IT IS ABSOLUTELY NECESSARY TO PERFORM THE FOLLOWING OPERATIONS:

- turn the appliance upside down and locate the small square door on the lower part;


- temporarily remove the flap;

Screw that must be unscrewed and removed

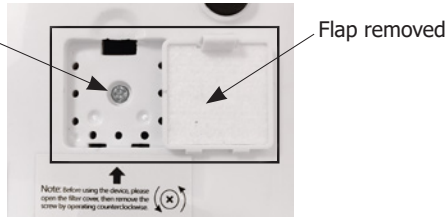- in the small compartment just opened, locate the screw that must be unscrewed and removed using the appropriate screwdriver supplied;

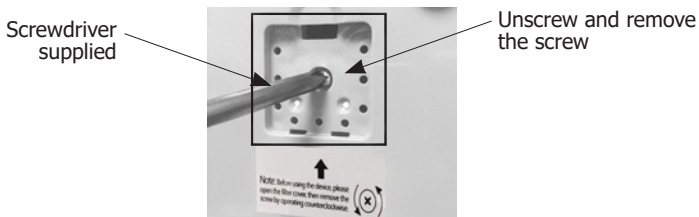
- once the screw has been removed, reposition the flap;
- recoil the equipment again by resting it on its base;
- proceed as described in the following paragraphs of the current chapter and in the chapter "Instructions for use of 4-chamber applicators."

Pressure

The pressure can be adjusted from 30 mmHg to 240 mmHg (+/- 10%)

To decrease the pressure,
turn the knob in a counter-clockwise direction.

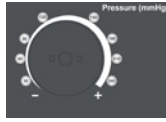To increase the pressure,
turn the knob in a clockwise direction.

- Make 2 inflations cycles to a minimum and then regulate the pressure.
- The inflation of each cuff takes approximately 4 seconds and then it switches to the next cuff.

Starting the massage at a minimum level of pressure. Gradually increase the pressure according to your own needs.

Time

The time can be adjusted from 0 to 30 minutes. After the unit has operated for a period of 30 minutes, it is necessary to wait at least 30 minutes (with the unit switched off) before starting a new session of operation. **Attention:** if you turn the knob on the "Always On" command, the automatic switch-off is deactivated; therefore the machine must be switched off manually.


To decrease the time,
turn the knob in a counter-clockwise direction.

To increase the times
turn the knob in a clockwise direction.

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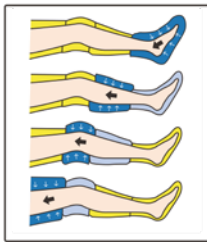Disposal as WEEE

Use exclusively inside

Electrical isolation class: Class II

Applicators and accessories

IMPORTANT: the applicators and the accessories provided with the machine may vary depending on the purchased configuration.

Leg cuffs

Leg cuff 4 chambers

The leg cuff/s are the accessories (shaped as a "leg") that must be worn in order to perform the massage on the lower limbs. They consist of **4 air chambers**.

In order to prevent knee pain, avoid bending your leg during the massage.

For user instructions, please refer to the following chapters.

Abdominal/buttocks belt

The abdominal belt is the accessory that must be worn in order to perform the massage on the abdomen and on the buttocks. It consists of **4 Standard air chambers**.

Do not place the belt on your breasts.

For user instructions, please refer to the following chapter.

Abdominal/buttocks belt
4 chambers

Arm cuff

Arm cuff 4 chambers

The arm cuff is the appropriate accessory composed of **4 air chambers** to wear for massaging the upper extremities.

In order to avoid that it slips down, secure the arm cuff using the appropriate Velcro brace.

For the instructions for use, please refer to the following sections.

Single connector

The single connector is an accessory made up of rubber tubes joined together that ends, at one end, with an attachment to be connected to the machine or to the splitter and, at the other, with some nozzles (air outlets) that must be inserted into the seats located on the applicator.

The single connector allows the use of the Mesis® PressoMassaggio® Plus + with 1 arm cuff or 1 leg cuff or the abdominal/buttock band or the 2 leg cuffs. In the latter case the configuration includes 2 single connectors plus a splitter.

For user instructions, please refer to the following chapters.


Single connector

Splitter

The splitter is an accessory made up of rubber tubes joined together with two ends of which one is a "male" connection with pins and the other is a twin "female" connection with holes. The splitter is combined with configurations with two or three applicators.

1. In the configuration with two applicators is contained 1 splitter. It allows two applicators to be used simultaneously such as 2 leg cuffs;

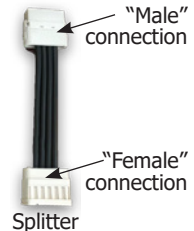
Splitter

2. In the configuration with three applicators there are 2 splitters. They allow three applicators to be used simultaneously such as 2 leg cuffs and the abdominal/buttocks belt.


In **case 1.** the "male" connection of the splitter must be connected to the machine and the "female" twin one, to the two single connectors of the two applicators.

In **case 2.** the two splitters will be connected to each other. Specifically, the 1st splitter will be connected to the machine as in case 1. While the "male" connection of the 2nd splitter will be inserted in one of the two seats of the "female" connection of the 1st splitter. A single connector must be connected to the second free seat of the "female" connection of the 1st splitter. In each of the two seats of the "female" connection of the 2nd splitter the remaining two single connectors will be connected.

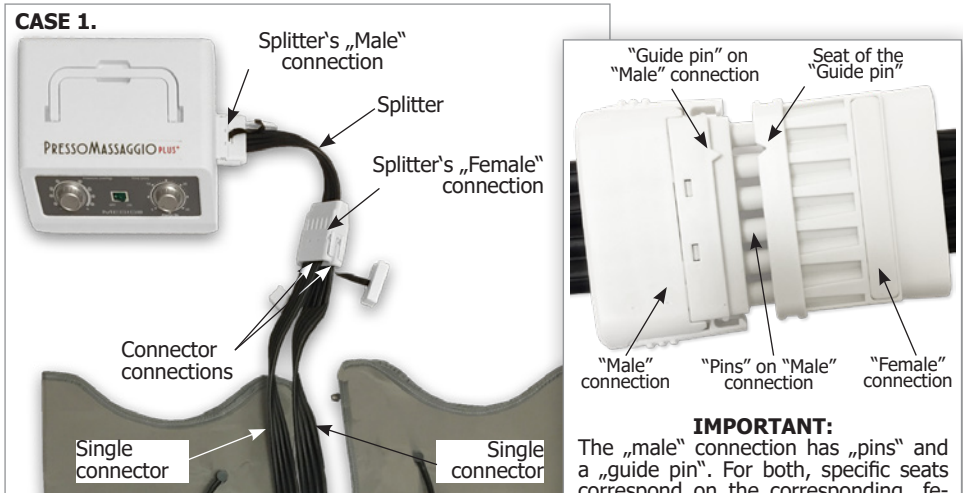IMPORTANT:

The „male“ connection has „pins“ and a „guide pin“. For both, specific seats correspond on the corresponding „female“ connection. Therefore pay close attention to the direction in which the „male“ and „female“ connections are connected together, since the „pins“ and „guide pin“ must be inserted in the corresponding seats.

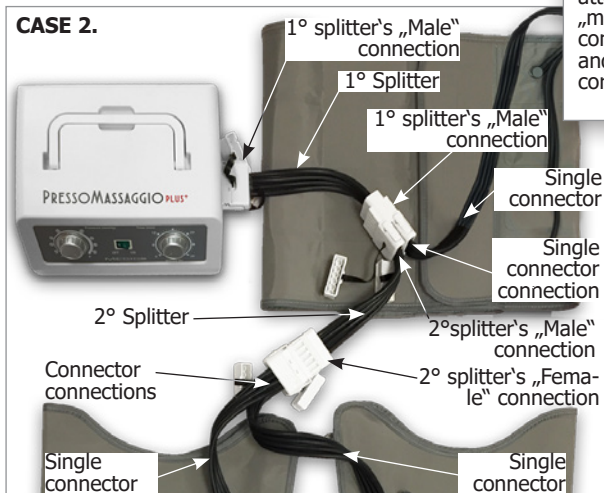Instructions for use applicators with 4 chambers

IMPORTANT NOTE: the applicators and the accessories provided with the machine may vary depending on the purchased model.

1. Connect the connector to the applicator that you want to use.

Choose the connector and/or the splitter to be connected depending on the applicators you would like to use:

- use the single connector with the arm cuff (Fig. A);
- use the single connector with the abdominal/buttocks belt (Fig. B);
- use the single connector with the leg cuff (Fig. C);
- use two single connectors and one splitter with two leg cuffs (Fig. D);
- use three single connectors and two splitter with two leg cuffs and the abdominal/buttocks belt (Fig. E);

Connect each nozzle (air outlet) on the end of the rubber hoses of the single connector to the relative inset on the applicator/s you want to use.

Nozzle (air outlet) on the end of the rubber hoses of connectors

Insert present on the applicators

ATTENTION: make the connection by combining the nozzle and the insert, according to the color. Nozzle and insert must have the same color.

Arm cuff, Abdominal/buttocks belt and Leg cuff used individually

- connect the nozzles on the end of the rubber tubes of the single connector to the insert on the applicator you want to use.

Fig. A

Fig. B

Fig. C

Leg cuffs

- if you want to use the two leg cuffs, connect the nozzles of the two single connectors to the insert on the legs. After that, connect the individual connectors to the "female" connection of the splitter.

ATTENTION: more detailed instructions on the connection between the connectors and the splitter on page 45 "Case 1."

Fig. D

Abdominal/buttocks belt

- if you wish to use the abdominal/buttocks belt and leg cuffs, connect the nozzles of the single connectors to the inserts on the belt and on the leg cuffs. After that connect the individual connectors of the leg cuffs to the "female" connection of the 2nd splitter and the "male" connection of the latter to the 1st splitter. Connect the single connector of the belt to the "female" connection of the 1st splitter.

ATTENTION: more detailed instructions on the connection between the connectors and the splitter on page 45 "Case 2."



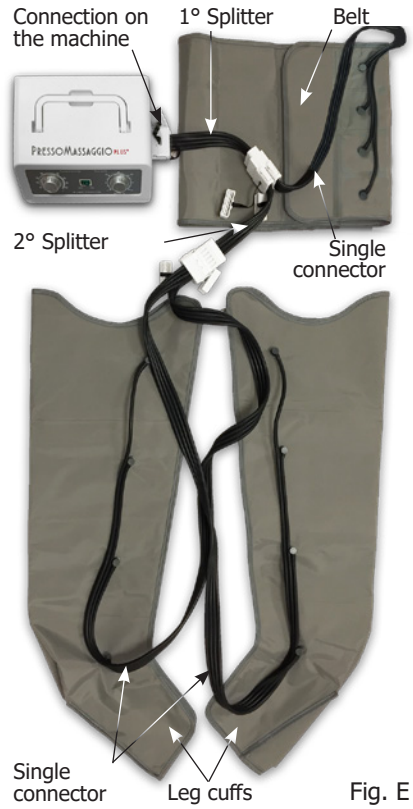
Fig. E

2. Connect the power cord.

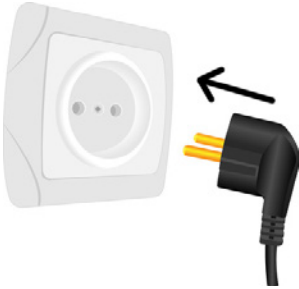3. Connect the single connector or the splitter, to the unit.

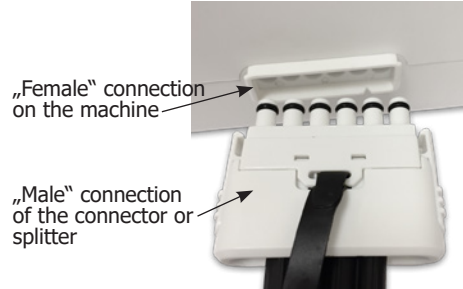In the PressoMassaggio® Plus⁺, the connection "machine-connector or divisor" is located on the side of the machine.

4. Put on the applicators that you wish to use.

Leg cuffs

- Put on the leg cuffs and close them with the zipper, place the leggings with the attachments turned inwards.
- **Inflation must start at the feet and gradually reach the groin.**


Abdominal/buttocks belt

- Put on abdominal belt on the buttocks-hips-abdomen and fasten it with the Velcro.
- **The inflation must start from the upper part of the abdomen and proceed downwards (groin).**

Arm cuff

- Put on the arm cuff.
- Fasten the velcro "strap" on the shoulder opposite the one with the arm cuff;
- **The inflation must start from the hand part and then go toward the shoulder in sequence.**

5. Press the power button to start the unit.

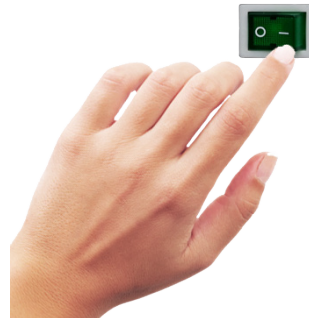In the PressoMassaggio® Plus⁺, the power button is located on the front of the machine.

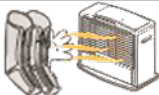
IMPORTANT: for the instructions for use of the unit, please refer to the relevant section.

IMPORTANT:

Make sure that the inflation takes place in the direction indicated on page 39 in the "Inflation direction" section of the "Features Mesis® PressoMassaggio® Plus⁺" chapter.

Maintenance and troubleshooting

The accessories must not be repaired and handled with care.

Do not position the leg cuffs, arm cuff, abdominal belt or other accessories near direct or indirect sources of heat.		
Do not place the device in damp or wet areas. Do not soak the equipment in any way, do not pour any liquid over the equipment.		
Do not wash the equipment and accessories with water or detergents. Simply wipe with a damp cloth.		
Do not store the device with its accessories in a particularly cold area or under a temperature of 0° C. The device could seriously be damaged.		
If you don't use the device for a length of time, place it back in its box.		

Environment of use and safekeeping

- Use and maintain the device and its accessories in a clean environment and at a temperature between 10 – 40 °C, and humidity of 30 – 85%;
- Do not use in outside environments.

Problem encountered	Possible solution
The unit does not turn on	<ul style="list-style-type: none"> • Check that it is correctly connected to the mains; • Check that the switch has been switch to ON (I) (OFF/ON (0-I) button); • Check the status of the power cable, it must be intact.
You hear strange hits similar to that of escaping air	<ul style="list-style-type: none"> • Check to ensure that the hoses, the accessories and the connector hoses are intact; • Check to ensure that all the connections (do the unit and to the applicators) have been inserted correctly; • Check that the air connector hoses are not bent, squashed or crushed.
The air does not flow to the applicators, there is insufficient pressure	<ul style="list-style-type: none"> • Check to ensure that the connector hoses to the unit and the applicators have been correctly connected; • Check that the air connector hoses are not bent, squashed or crushed; • Check that the chambers are inflated progressively and according to the selected program; • Make sure you set the pressure control correctly; • Verify that the actual air outlet from the unit is present.
The air continues to flow into the chambers and you hear an air leak in the applicators	<ul style="list-style-type: none"> • Check for holes or damage in the applicators' air chambers. If no abnormalities are present, switch off the device, turn it on again to check again to ensure that is functioning properly.

Warranty

- The device is guaranteed for the first user for a period of twenty four (24) months, three (3) months for the applicators and accessories subject to wear from the date of purchase against defects in materials or manufacture, provided that it is used properly and maintained in good working condition.

The validity of the warranty shall be limited in the following cases:

- There is no warranty for accessories and materials that are considered to be „consumables“.
- The present warranty shall be valid and effective in the country where the product has been purchased. In the event the product is purchased in any country belong to the European Union, the warranty shall then be valid in all of its member states.
- Problems due to having dropped the device or improper voltage use.
- If the device has been opened or modified by unauthorised personnel.
- If the device has been used in unsuitable environments.
- If the device has not been used in accordance with the instructions contained in this manual.
- Warranty on the device: 2 years.
- Warranty for the accessories and applicators: 3 months.

To obtain the warranty service, the user must observe the following warranty clauses:

- The products must be delivered for repair in the original packaging.
- The product warranty will be subject to the production of a tax document (sales receipt, payment receipt, sales invoice or mail-order sales document via an online shop or catalogue), attesting the product's date of purchase.
- Repair will not be performed on the original date of the warranty's expiry and it will not constitute any renewal or extension of the same.
- If, at the time of repair, there is no fault in the product found, expenses related to the technical verification shall be charged.
- The warranty shall lapse if malfunction has been determined to be due to impact, falls, incorrect or improper use of the product, use of a power supply or charger that is not the original, accidents, alterations to the product, substitution/removal of the warranty seals and/or tampering with the product.
- Moreover, the same warranty shall not cover damage caused during transport through use of unsuitable packaging.
- The warranty does not cover the impossibility to use the product, other incidental or consequential costs or other expenses incurred by the purchaser.

Product support

In the case of breakdown or malfunction of the device, contact your dealer.

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Prodotto da: Dili (Xiamen) Import And Export Co., Ltd